


Assay Class: High Sensitivity DNA Assay  
Data Path: T:\...gh Sensitivity DNA Assay\_DE72902247\_2020-04-20\_12-02-44.xad

Created: 4/20/2020 12:02:44 PM  
Modified: 4/20/2020 12:39:01 PM

### Electrophoresis File Run Summary

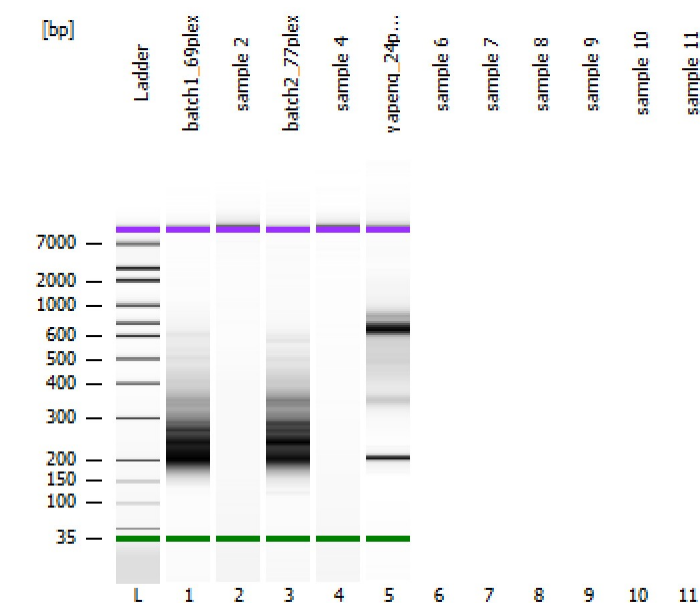

#### Instrument Information:


Instrument Name: DE72902247      Firmware: C.01.069  
Serial#: DE72902247      Type: G2939A

#### Assay Information:

Assay Origin Path: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy  
Assay Class: High Sensitivity DNA Assay  
Version: 1.03  
Assay Comments: Copyright © 2003-2010 Agilent Technologies

#### Chip Information:

Chip Lot #:  
Reagent Kit Lot #:  
Chip Comments:

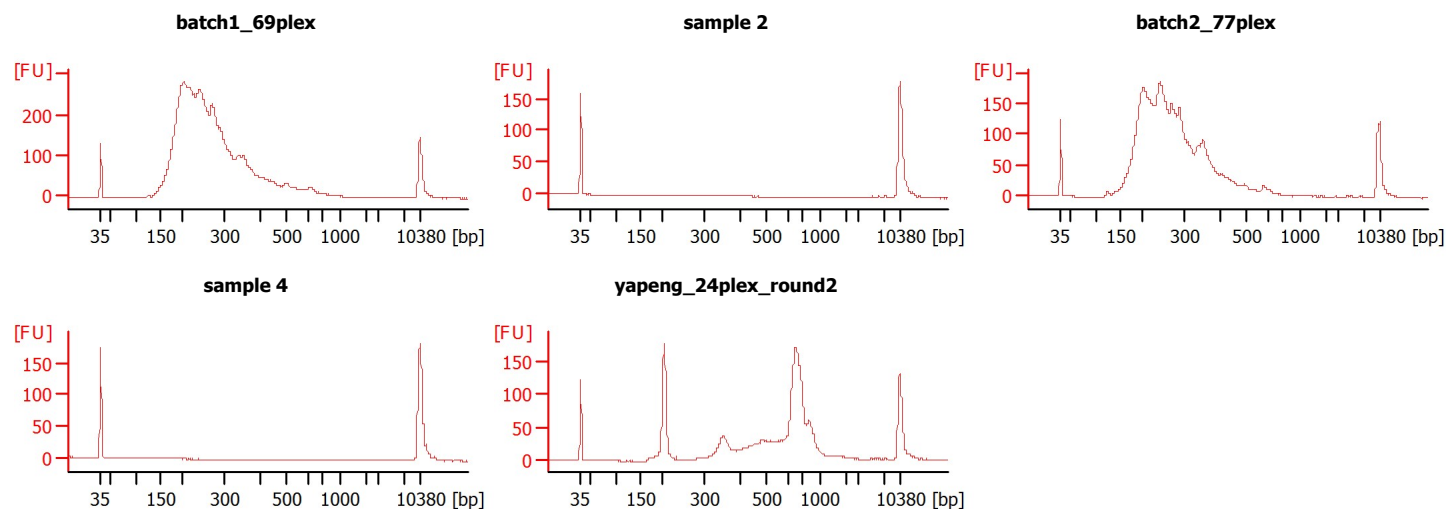

Assay Class: High Sensitivity DNA Assay  
Data Path: T:\...gh Sensitivity DNA Assay\_DE72902247\_2020-04-20\_12-02-44.xad

Created: 4/20/2020 12:02:44 PM  
Modified: 4/20/2020 12:39:01 PM

**Electrophoresis File Run Summary (Chip Summary)**

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
batch1_69plex		<input type="checkbox"/>	✓			
sample 2		<input type="checkbox"/>	✓			
batch2_77plex		<input type="checkbox"/>	✓			
sample 4		<input type="checkbox"/>	✓			
yapeng_24plex_round 2		<input type="checkbox"/>	✓			
sample 6		<input type="checkbox"/>				
sample 7		<input type="checkbox"/>				
sample 8		<input type="checkbox"/>				
sample 9		<input type="checkbox"/>				
sample 10		<input type="checkbox"/>				
sample 11		<input type="checkbox"/>				
Ladder		<input type="checkbox"/>	✓			

**Chip Lot #****Reagent Kit Lot #****Chip Comments :**

Assay Class: High Sensitivity DNA Assay  
Data Path: T:\...gh Sensitivity DNA Assay\_DE72902247\_2020-04-20\_12-02-44.xad

Created: 4/20/2020 12:02:44 PM  
Modified: 4/20/2020 12:39:01 PM

## Electrophoresis Assay Details

### General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12  
Minimum Visible Range [s] : 32  
Maximum Visible Range [s] : 138  
Start Analysis Time Range [s] : 33  
End Analysis Time Range [s] : 137.5  
Ladder Concentration [pg/ $\mu$ l] : 1950  
Uses Standard Area for Ladder Fragments  
Lower Marker Concentration [pg/ $\mu$ l] : 125  
Upper Marker Concentration [pg/ $\mu$ l] : 75  
Used Upper Marker for Quantitation  
Standard Curve Fit is Point to Point  
Show Data Aligned to Lower and Upper Marker

### Integrator Settings

Integration Start Time [s] : 33.05  
Integration End Time [s] : 137  
Slope Threshold : 0.8  
Height Threshold [FU] : 5  
Area Threshold : 0.1  
Width Threshold [s] : 0.6  
Baseline Plateau [s] : 0.5

### Filter Settings

Filter Width [s] : 0.5  
Polynomial Order : 4


### Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay  
 Data Path: T:\...gh Sensitivity DNA Assay\_DE72902247\_2020-04-20\_12-02-44.xad

Created: 4/20/2020 12:02:44 PM  
 Modified: 4/20/2020 12:39:01 PM

### Electropherogram Summary

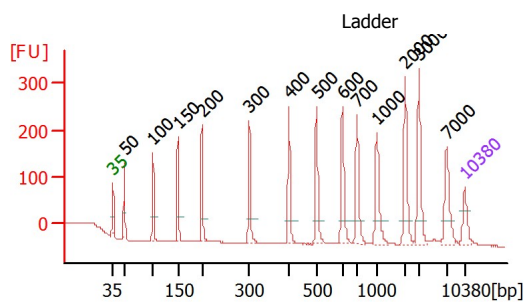

#### Overall Results for Ladder

Noise: 0.3


#### Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	50	150.00	4,545.5	Ladder Peak
3	100	150.00	2,272.7	Ladder Peak
4	150	150.00	1,515.2	Ladder Peak
5	200	150.00	1,136.4	Ladder Peak
6	300	150.00	757.6	Ladder Peak
7	400	150.00	568.2	Ladder Peak
8	500	150.00	454.5	Ladder Peak
9	600	150.00	378.8	Ladder Peak
10	700	150.00	324.7	Ladder Peak
11	1,000	150.00	227.3	Ladder Peak
12	2,000	150.00	113.6	Ladder Peak
13	3,000	150.00	75.8	Ladder Peak

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
14	7,000	150.00	32.5	Ladder Peak
15	10,380	75.00	10.9	Upper Marker

Assay Class: High Sensitivity DNA Assay  
Data Path: T:\...gh Sensitivity DNA Assay\_DE72902247\_2020-04-20\_12-02-44.xad

Created: 4/20/2020 12:02:44 PM  
Modified: 4/20/2020 12:39:01 PM

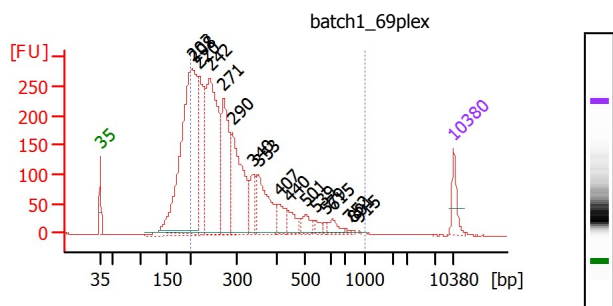**Electropherogram Summary Continued ...****Overall Results for sample 1 : batch1\_69plex**


Number of peaks found: 17  
Noise: 0.4  
Corr. Area 1: 5,221.0

**Region table for sample 1 : batch1\_69plex**

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarit y [pmol/l]	Co lo r
200	1,000	5,221.0	81	301	34.9	4,348.33	24,333.4	■

**Peak table for sample 1 :**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	207	6,389.13	46,766.0	
3	208	1,745.85	12,695.5	
4	220	570.39	3,925.5	
5	242	1,072.53	6,704.3	
6	271	569.74	3,191.0	
7	290	569.03	2,977.3	
8	340	155.13	691.7	
9	353	329.00	1,410.4	
10	407	99.76	371.7	
11	440	73.76	253.9	
12	501	67.05	202.7	
13	539	31.16	87.6	
14	579	16.37	42.8	
15	615	47.51	117.1	
16	752	4.41	8.9	
17	801	3.88	7.3	
18	915	2.04	3.4	
19	10,380	75.00	10.9	Upper Marker

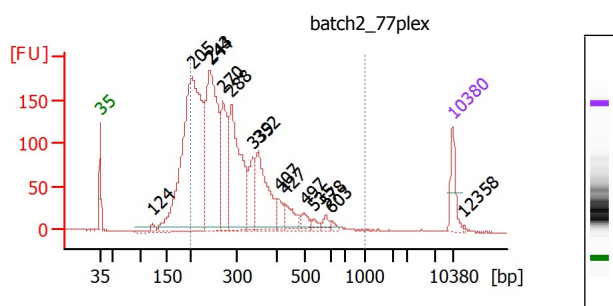**Overall Results for sample 3 : batch2\_77plex**

Number of peaks found: 15  
Noise: 0.3  
Corr. Area 1: 3,535.2

**Region table for sample 3 : batch2\_77plex**

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarit y [pmol/l]	Co lo r
200	1,000	3,535.2	83	299	31.1	3,140.93	17,416.6	■

**Peak table for sample 3 :**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	243	3,996.45	24,933.3	
3	124	17.72	216.6	
4	205	1,413.70	10,473.3	
5	244	836.17	5,187.2	
6	270	348.50	1,958.2	
7	288	522.17	2,750.3	
8	339	139.67	624.2	
9	352	307.54	1,324.7	
10	407	55.98	208.3	
11	427	73.19	259.7	
12	497	36.08	110.0	
13	535	17.74	50.3	
14	578	20.70	54.2	
15	603	11.24	28.2	
16	10,380	75.00	10.9	Upper Marker
17	12,358	0.00	0.0	

Assay Class: High Sensitivity DNA Assay  
Data Path: T:\...gh Sensitivity DNA Assay\_DE72902247\_2020-04-20\_12-02-44.xad

Created: 4/20/2020 12:02:44 PM  
Modified: 4/20/2020 12:39:01 PM

**Gel Image**